# **Elabscience**®

## Elab Fluor<sup>®</sup> 488 Anti-Mouse/Rat CD29 Antibody[HMβ1-1]

Catalog No.E-AB-F1309LStorageStore at 2~8°C, Avoid freeze / thaw cycles

Reactivity Applications

Mouse,Rat FCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

#### **Antigen Information**

Alternate Names	integrin β1, VLA-β chain, β1 integrin, GPIIa, ITGB1
Uniprot ID	P09055, P49134
Gene ID	16412, 24511
Background	CD29 is a 130 kD protein, also known as integrin $\beta$ 1, VLA- $\beta$ chain, or GPIIa. It is a member of the integrin family, expressed broadly on leukocytes, endothelial cells, smooth muscle, and epithelial cells. In association with CD49a-f, CD29 forms the VLA-1 through VLA-6 complexes, respectively. It plays an important role in cell-cell or cell-matrix interaction. The HM $\beta$ 1-1 antibody reacts with both mouse and rat CD29. It is able to block cell adhesion and inhibit T cell proliferation.

#### **Product Details**

Form	Liquid
Size	50Tests/100Tests/100Tests×2
Clone No.	ΗΜβ1-1
Host	Armenian Hamster
Isotype	Armenian Hamster IgG
Reactivity	Mouse,Rat
Application	FCM
Isotype Control	[Product E-AB-F09852L]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

**For Research Use Only** 

# **Elabscience**®

### Fluorophore


#### **Conjugation:** Elab Fluor<sup>®</sup> 488

Elab Fluor<sup>®</sup> 488 is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 520 nm (e.g., a 525/40 nm bandpass filter).

#### **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

### **Product data**

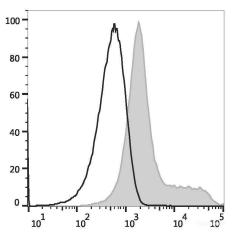

C57BL/6 murine splenocytes are stained with Elab Fluor<sup>®</sup> 488 Anti-Mouse/Rat CD29 Antibody (filled gray histogram) or Elab Fluor<sup>®</sup> 488 Armenian Hamster IgG Isotype Control (empty black histogram).

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry <u>https://www.elabscience.com/List-detail-5594.html</u>
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>