



A Reliable Research Partner in Life Science and Medicine

PE Anti-Human CD79a Antibody[HM47]

Catalog No. E-AB-F1370D Reactivity Human Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications ICFCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Mb-1, Iga, CD79

Uniprot ID P11912 Gene ID 973

Background CD79a is a 47 kD type I integral membrane protein, also known as mb-1 or Iga. It is a member of

the Ig superfamily and disulphide-associated with CD79b (B29). The interaction of

CD79a/CD79b heterodimer with B cell suface Ig forms B cell antigen receptor complex. CD79a is expressed in B cells from early pre-B to plasma cell stage. It has been shown that CD79a is also weakly expressed in some precursors of T- and myeloid cells. CD79 mediates the transport of

IgM to B cell surface and transduces signals initiated by BCR aggregation.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. **HM47** Host Mouse

Mouse IgG1, κ **Isotype**

Reactivity Human **Application ICFCM**

Isotype Control PE Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792D]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017 Web: www.elabscience.com

Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: PE

PE is designed to be excited by the Blue (488 nm), Green (532 nm) and Yellow-Green (561 nm) lasers and detected using an optical filter centered near 575 nm (e.g., a 585/42 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

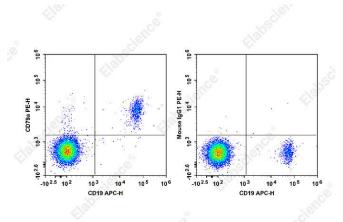

Staining of normal human peripheral blood cells with APC Anti-Human CD19 Antibody and PE Anti-Human CD79a Antibody[HM47] (left) or PE Mouse IgG1, κ Isotype Control (right). Cells in the lymphocytes gate were used for analysis.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com