

Recombinant Human TIMP2/TIMP-2 Protein (His Tag)

Catalog No. PKSH033118

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Synonyms	Metalloproteinase Inhibitor 2;CSC-21K;Tissue Inhibitor of Metalloproteinases 2;CSC-21K;DDC8
Species	Human
Expression Host	HEK293 Cells
Sequence	Cys27-Pro220
Accession	P16035
Calculated Molecular Weight	22.8 kDa
Observed molecular weight	26 kDa
Tag	C-His
Bioactivity	Not validated for activity

Properties

Purity	> 95 % as determined by reducing SDS-PAGE.
Endotoxin	< 1.0 EU per µg of the protein as determined by the LAL method.
Storage	Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of reconstituted samples are stable at < -20°C for 3 months.
Shipping	This product is provided as lyophilized powder which is shipped with ice packs.
Formulation	Lyophilized from a 0.2 µm filtered solution of 20mM PB, 150mM NaCl, pH 7.2. Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as protectants before lyophilization. Please refer to the specific buffer information in the printed manual.
Reconstitution	Please refer to the printed manual for detailed information.

Data

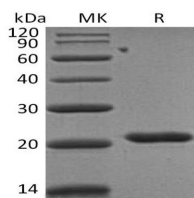

> 95 % as determined by reducing SDS-PAGE.

Background

Tissue inhibitors of metalloproteinases or TIMPs are a family of proteins that regulate the activation and proteolytic

For Research Use Only

activity of the zinc enzymes known as matrix metalloproteinases (MMPs). There are four members of the family, TIMP-1, TIMP-2, TIMP-3, and TIMP-4. Tissue Inhibitor of Metalloproteinases 2 (TIMP-2) is a non N-glycosylated protein with a molecular mass of 22 kDa. It produced by a wide range of cell types, which inhibits MMPs non-covalently by the formation of binary complexes and irreversibly inactivates them by binding to their catalytic zinc cofactor. TIMP-2 also has erythroid-potentiating and cell growth promoting activities.