Elabscience®

FITC Anti-Human CD5 Antibody[UCHT2]

Catalog No. E-AB-F1041C Storage Store at 2~8°C, Avoid freeze / thaw cycles Reactivity Human Applications

FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

| Alternate Names | T-cell surface glycoprotein CD5, CD5, Lymphocyte antigen T1/Leu-1, LEU1 |
|-----------------|---|
| Uniprot ID | P06127 |
| Background | CD5 is a 67 kD single chain type I glycoprotein also known as Leu-1, Ly-1 and T1. It is a member of the scavenger receptor superfamily found on T cells, thymocytes, B cell subsets, chronic B lymphocytic leukemia (B-Cells), and peripheral blood dendritic cells. CD5 modulates T and B cell receptor signaling, thymocyte maturation, and T-B cell interactions upon binding to ligands such as CD72. |

Product Details

| Form | Liquid |
|---------------------|---|
| Size | 20Tests/100Tests/100Tests×2 |
| Clone No. | UCHT2 |
| Host | Mouse |
| Isotype | Mouse IgG1, κ |
| Reactivity | Human |
| Application | FCM |
| Isotype Control | FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C] |
| Storage Buffer | Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. |
| Shipping | Biological ice pack at 4 °C |
| Stability & Storage | Keep as concentrated solution. |
| | Store at 2~8°C and protected from prolonged exposure to light.Do not freeze. |
| | This product is guaranteed up to one year from purchase. |

For Research Use Only

Elabscience®

Fluorophore


Conjugation: FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

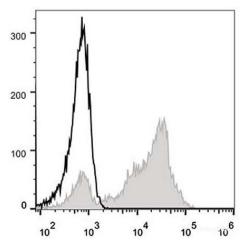

Human peripheral blood lymphocytes are stained with FITC Anti-Human CD5 Antibody (filled gray histogram). Unstained lymphocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry <u>https://www.elabscience.com/List-detail-5594.html</u>
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>