

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 647 Anti-Mouse CD19 Antibody[1D3]

E-AB-F0986M Catalog No. Reactivity Mouse Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names B-lymphocyte antigen CD19, Cd19, Differentiation antigen CD19, CD19

Uniprot ID

Background CD19 is a 95 kD glycoprotein also known as B4. It is a member of the Ig superfamily, expressed

> on all pro-B to mature B cells (during development) and follicular dendritic cells. Plasma cells do not express CD19. CD19, in association with CD21 and CD81, forms a molecular complex

integral to B cell activation.

Product Details

Form Liquid

50Tests/100Tests/100Tests×2 Size

Clone No. 1D3 Host Rat

Isotype Rat IgG2a, κ Mouse Reactivity **Application FCM**

Isotype Control Elab Fluor[®] 647 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09832M]

Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. **Storage Buffer**

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Email: techsupport@elabscience.com Web: www.elabscience.com

Tel: 1-832-243-6086

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® 647

Elab Fluor $^{\odot}$ 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

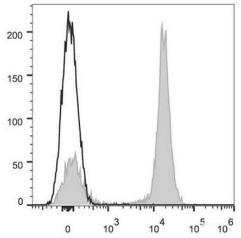

C57BL/6 murine splenocytes are stained with Elab Fluor[®] 647 Anti-Mouse CD19 Antibody (filled gray histogram). Unstained splenocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com