

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 488 Anti-Mouse CD86 Antibody[GL-1]

E-AB-F0994L Catalog No. Reactivity Mouse Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names T-lymphocyte activation antigen CD86,Cd86,Activation B7-2 antigen,Early T-cell costimulatory

molecule 1,ETC-1

Uniprot ID P42082

Background CD86 is an 80 kD immunoglobulin superfamily member also known as B7-2, B70, and Ly-58.

> CD86 is expressed on activated B and T cells, macrophages, dendritic cells, and astrocytes. CD86, along with CD80, is a ligand of CD28 and CD152 (CTLA-4). CD86 is expressed earlier in the immune response than CD80. CD86 has also been shown to be involved in immunoglobulin classswitching and triggering of NK cell-mediated cytotoxicity. CD86 binds to CD28 to transduce costimulatory signals for T cell activation, proliferation, and cytokine production. CD86 can also

bind to CD152, also known as CTLA-4, to deliver an inhibitory signal to T cells.

Product Details

Form Liquid

50Tests/100Tests/100Tests×2 Size

Clone No. GL-1 Host Rat

Rat IgG2a, κ **Isotype** Reactivity Mouse **Application FCM**

Elab Fluor[®] 488 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09832L] **Isotype Control**

Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. **Storage Buffer**

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Email: techsupport@elabscience.com Web: www.elabscience.com



A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® 488

Elab Fluor[®] 488 is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 520 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

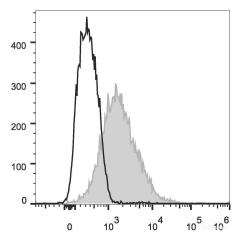

C57BL/6 murine splenocytes are stained with Elab Fluor[®] 488 Anti-Mouse CD86 Antibody (filled gray histogram). Unstained splenocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com