



A Reliable Research Partner in Life Science and Medicine

APC Anti-Human CD58 Antibody[TS2/9.1]

E-AB-F1068E Catalog No. Reactivity Human Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Uniprot ID

Lymphocyte function-associated antigen 3,CD58,Ag3,Surface glycoprotein LFA-3,LFA3

Background CD58, also known as lymphocyte function-associated antigen 3 (LFA-3) is a 45-70 kD cell

surface protein that is a member of the immunoglobulin superfamily. Alternative splicing of CD58 gives rise to transmembrane and glycosylphosphatidylinositol (GPI)-anchored forms on cell surface. CD58 is expressed on both hematopoietic and non-hematopoietic cells including B cells, T cells, monocytes, erythrocytes, endothelial cells, epithelial cells, and fibroblasts. High levels are observed on memory T cells and dendritic cells. CD58 expressed on antigen presenting cells and target cells enhances T cell recognition via the binding of it's cognate ligand, CD2, on the T cell

surface.

Product Details

Form Liquid

20Tests/100Tests/100Tests×2 Size

Clone No. TS2/9.1 Host Mouse Mouse IgG1, κ **Isotype**

Reactivity Human **Application FCM**

Isotype Control APC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792E]

Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. **Storage Buffer**

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017 Web: www.elabscience.com

Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: APC

APC is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 660 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

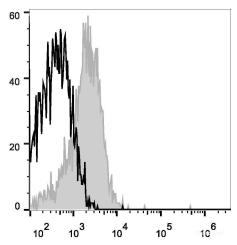

Human peripheral blood lymphocytes are stained with APC Anti-Human CD58 Antibody (filled gray histogram).

Unstained lymphocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com