



A Reliable Research Partner in Life Science and Medicine

APC Anti-Human CD10 Antibody[CB-CALLA]

Catalog No.E-AB-F1078EReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Neprilysin,MME,CALLA,NEP,Neutral endopeptidase,SFE,CD10

Uniprot ID P08473

Background CD10 is a 100 kD neutral endopeptidase and a member of the metalloprotease family. It is a type

II transmembrane protein also known as common acute lymphoblastic leukemia antigen (CALLA), enkephalinase, and neprilysin. CD10 is expressed on B cell precursors, T cell

precursors, and neutrophils. CD10 is involved in B cell development and has been shown to bind opioid enkephalins, bradykinin, angiotensins I and II, and other biologically active peptides.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. CB-CALLA Host Mouse

Isotype Mouse IgG1, κ

Reactivity Human **Application** FCM

Isotype Control APC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792E]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com

Tel: 1-832-243-6086 Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: APC

APC is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 660 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

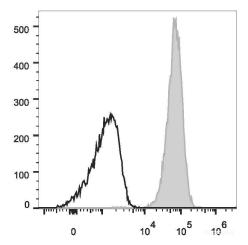

Human peripheral blood lymphocytes are stained with APC Anti-Human CD10 Antibody (filled gray histogram).

Unstained lymphocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com