



A Reliable Research Partner in Life Science and Medicine

FITC Anti-Human CD9 Antibody[HI9a]

E-AB-F1086C Catalog No. Reactivity Human Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names CD9 antigen, CD9,5H9 antigen, Cell growth-inhibiting gene 2 protein, Leukocyte antigen

MIC3,MRP-1,Tspan-29,p24,CD9,MIC3, TSPAN29

Uniprot ID P21926

Background CD9 is a 24 kD type III transmembrane protein also known as tetraspanin, MRP-1 and DRAP-24.

It is a member of the tetraspan family (spanning the membrane four times) found on platelets, B cell progenitors, activated lymphocytes, granulocytes, endothelial cells and epithelial cells. CD9 induces adhesion, platelet aggregation, and B cell development. CD9 has been shown to associate

with CD63, CD81, CD82, and CD36 and to bind to β1 integrins.

Product Details

Form Liquid

20Tests/100Tests/100Tests×2 Size

Clone No. HI9a Mouse Host

Isotype Mouse IgG1, κ

Human Reactivity **Application FCM**

Isotype Control FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Biological ice pack at 4 °C **Shipping** Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Email: techsupport@elabscience.com Web: www.elabscience.com

Tel: 1-832-243-6086 Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

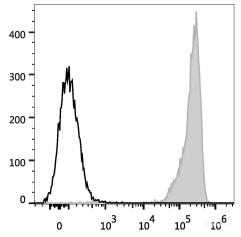

Human platelets are stained with FITC Anti-Human CD9 Antibody (filled gray histogram). Unstained platelets (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com