# **Elabscience**®

### PE/Cyanine5 Anti-Human CD41 Antibody[HIP8]

Catalog No.E-AB-F1088GStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityApplications

Human FCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

### **Antigen Information**

Alternate Names	Integrin alpha-Iib,ITGA2B,GPalpha Iib,GPIIb,Platelet membrane glycoprotein Iib,CD41,GP2B, ITGAB
Uniprot ID	P08514
Background	CD41 is a 125/25 kD $\alpha$ subunit of the gpIIb/IIIa (CD41/CD61) complex. CD41 is a heterodimer composed of a heavy chain (gpIIb $\alpha$ ) and light chain (gpIIb $\beta$ ) linked by a single disulfide bond. It is a member of the integrin family primarily expressed on platelets and megakaryocytes. CD41 has been reported to be involved with platelet aggregation and platelet attachment to the ECM. CD41/CD61 complex acts as the receptor for fibrinogen, fibronectin, Von Willebrand factor and thrombin.

#### **Product Details**

Form	Liquid
-	1
Size	20Tests/100Tests/100Tests×2
Clone No.	HIP8
Host	Mouse
Isotype	Mouse IgG1, ĸ
Reactivity	Human
Application	FCM
Isotype Control	[Product E-AB-F09792G]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

**For Research Use Only** 

# **Elabscience**®

### Fluorophore


#### Conjugation: PE/Cyanine5

PE/Cyanine5 is designed to be excited by the Blue (488 nm), Green (532 nm) and yellow-green (561 nm) lasers and detected using an optical filter centered near 670 nm (e.g., a 690/50 nm bandpass filter).

#### **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

#### **Product data**

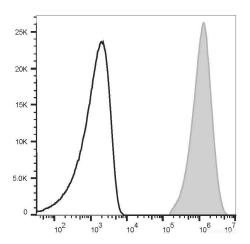

 $\label{eq:HIP8} Human \ platelets \ are \ stained \ with \ PE/Cyanine5 \ Anti-Human \ CD41 \ Antibody [HIP8] \ (filled \ gray \ histogram) \ or \ PE/Cyanine5 \ Mouse \ IgG1, \ \kappa \ Isotype \ Control \ (empty \ black \ histogram).$ 

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry <u>https://www.elabscience.com/List-detail-5594.html</u>
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>