

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

APC Anti-Mouse CD71 Antibody[R17 217.1.3/TIB-219]

Catalog No.E-AB-F1093EReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Transferrin receptor protein 1,Tfrc,TR,TfR,TfR1,Trfr,CD71

Uniprot ID Q6235

Background CD71 is a 95 kD type II heterodimeric transmembrane glycoprotein that is also known as T9 and

transferrin receptor. CD71 is expressed on proliferating cells, reticulocytes, and erythroid precursors. Its expression is very low on resting leukocytes. CD71 plays a role in the control of cellular proliferation by facilitating the uptake of iron via ferrotransferrin binding and the

recycling of apotransferrin to the cell surface.

Product Details

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. R17 217.1.3/TIB-219

Host Rat

 $\begin{array}{ll} \textbf{Isotype} & \text{Rat IgG2a, } \kappa \\ \textbf{Reactivity} & \text{Mouse} \\ \textbf{Application} & \text{FCM} \\ \end{array}$

Isotype Control APC Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09832E]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086

Web: www.elabscience.com
Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: APC

APC is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 660 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

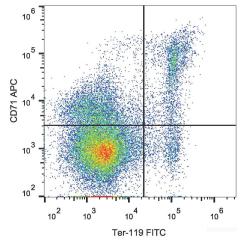

C57BL/6 mouse bone marrow cells were stained with FITC Anti-Mouse TER-119 Antibody and APC Anti-Mouse CD71 Antibody.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com