

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

PE/Cyanine7 Anti-Mouse CD49b/pan-NK cells Antibody[DX5]

Catalog No.E-AB-F1116HReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Integrin alpha-2,CD49 antigen-like family member B,Collagen receptor,Platelet membrane

glycoprotein Ia, GPIa, VLA-2 subunit alpha, CD49b

Uniprot ID Q62469

Background DX5 antigen has been recently characterized as CD49b. It is a 150 kD integrin α chain also

known as $\alpha 2$ integrin, VLA-2 α chain, and integrin $\alpha 2$ chain. CD49b non-covalently associates with CD29 ($\beta 1$ integrin) to form the CD49b/CD29 complex known as VLA-2, a receptor for collagen and laminin. CD49b is expressed on platelets, the majority of NK cells, NKT cells, and a small subset of CD8+ T cells (this population can be significantly increased following viral infection). DX5 is used for the identification and isolation of NK cells, and is especially useful

for identifying NK cells in mice lacking the NK1.1 antigen.

Product Details

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No.DX5HostRatIsotypeRat IgM, κReactivityMouseApplicationFCM

Isotype Control [Product E-AB-F09772H]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at $2\hbox{--}8^\circ C$ and protected from prolonged exposure to light. Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086

Web: www.elabscience.com
Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: PE/Cyanine7

PE/Cyanine7 is designed to be excited by the Blue (488 nm), Green (532 nm) and yellow-green (561 nm) lasers and detected using an optical filter centered near 775 nm (e.g., a 780/60 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

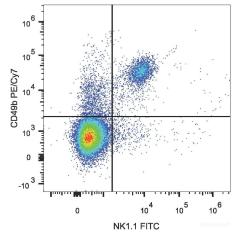

C57BL/6 murine splenocytes are stained with PE/Cyanine7 Anti-Mouse CD49b Antibody and FITC Anti-Mouse CD161/NK1.1 Antibody.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com