# **Elabscience**®

# FITC Anti-Human CD34 Antibody[581]

Catalog No.E-AB-F1143CStorageStore at 2~8°C, Avoid freeze / thaw cycles

Reactivity Applications

Human FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

### **Antigen Information**

Alternate Names Uniprot ID Background	Gp105-120, My10 P28906 CD34, also known as gp105-120, is a type I monomeric sialomucin-like glycophosphoprotein with an approximate molecular weight of 105-120 kD. Selectively expressed on the majority of hematopoietic stem/progenitor cells, bone marrow stromal cells, capillary endothelial cells, embryonic fibroblasts, and some nervous tissue, CD34 is a commonly used marker to identify human hematopoietic stem/progenitor cells. According to the differential sensitivity to enzymatic cleavage, four groups of epitopes of CD34 have been described. CD34 mediates cell adhesion and hymphoantes homing through hinding to L calactin ligende.
	lymphocytes homing through binding to L-selectin and E-selectin ligands.

#### **Product Details**

Form	Liquid
Size	20Tests/100Tests/100Tests×2
Clone No.	581
Host	Mouse
Isotype	Mouse IgG1, κ
Reactivity	Human
Application	FCM
Isotype Control	FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

**For Research Use Only** 

# **Elabscience**®

## Fluorophore


#### **Conjugation:** FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

#### **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

### **Product data**

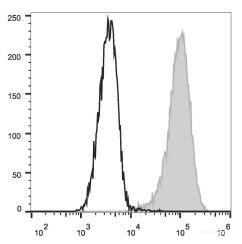

KG1 cells are stained with FITC Anti-Human CD34 Antibody (filled gray histogram). Unstained cells (empty black histogram) are used as control.

### **Related Information**

- 1. Sample Preparation for Flow Cytometry <u>https://www.elabscience.com/List-detail-5594.html</u>
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>