



A Reliable Research Partner in Life Science and Medicine

APC Anti-Human CD122/IL-2RB Antibody[TU27]

Catalog No.E-AB-F1151EReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Interleukin-2 receptor subunit beta,IL-2RB,p75,CD122,IL-2Rβ

Uniprot ID P14784

Background CD122 is a 70-75 kD type I transmembrane glycoprotein and member of the Ig superfamily. It is

IL-2 receptor β chain also known as IL-2R β , which is also shared by the IL-15 receptor. CD122 is constitutively expressed by NK cells and at lower levels by a subset of T cells. Its expression is upregulated upon activation. The IL-2R β chain can combine with either the common γ subunit (γ c, CD132) alone or with the γ c subunit and the IL-2R α subunit (CD25) to generate intermediate

or high affinity IL-2 receptor complexes, respectively. CD122 expression levels can be

upregulated by activation.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. TU27 Host Mouse

Isotype Mouse IgG1, κ

Reactivity Human **Application** FCM

Isotype Control APC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792E]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086
 Fax: 1-832-243-6017

 Web: www.elabscience.com
 Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: APC

APC is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 660 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

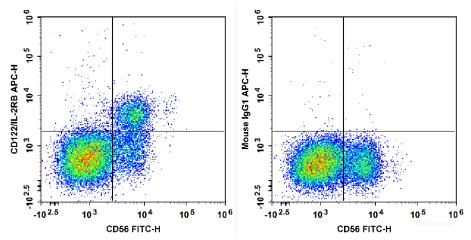

Human peripheral blood lymphocytes are stained with FITC Anti-Human CD56 Antibody and APC Anti-Human CD122/IL-2RB Antibody[TU27] (Left). Lymphocytes are stained with FITC Anti-Human CD56 Antibody and APC Mouse IgG1, κ Isotype Control (Right).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com