



A Reliable Research Partner in Life Science and Medicine

# APC Anti-Human CD36 Antibody[5-271]

Catalog No. E-AB-F1164E Reactivity Human Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM** 

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

## **Antigen Information**

**Alternate Names** GP3B,GP4,GPIV,GPIIIB,PAS-4

**Uniprot ID** P16671

**Background** CD36 is an 85 kD integral membrane glycoprotein, also known as GPIIIb, or GPIV. It is

> expressed on various epithelial and endothelial cells as well as erythrocytes, platelets, macrophages/monocytes and some macrophage-derived dendritic cells. CD36 functions as a scavenger receptor, binding thrombospondin, long chain fatty acids, oxidized LDL, collagen type I, IV, and V as well as apoptotic cells. The 5-271 antibody has been reported to be useful for flow

cytometry.

#### **Product Details**

Form Liquid

20Tests/100Tests/100Tests×2 Size

Clone No. 5-271 Mouse Host

**Isotype** Mouse IgG2a, κ

Human Reactivity **Application FCM** 

**Isotype Control** APC Mouse IgG2a, κ Isotype Control[C1.18.4] [Product E-AB-F09802E]

Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. **Storage Buffer** 

Biological ice pack at 4 °C Shipping Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017 Web: www.elabscience.com

Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

## **Fluorophore**


Conjugation: APC

APC is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 660 nm (e.g., a 660/20 nm bandpass filter).

## **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

### **Product data**

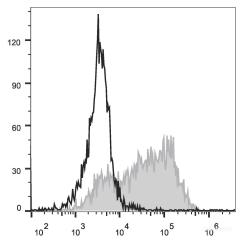

Human ErythroLeukemia cell line HEL are stained with APC Anti-Human CD36 Antibody (filled gray histogram). Unstained Human ErythroLeukemia cell line HEL (empty black histogram) are used as control.

### **Related Information**

- 1. Sample Preparation for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5594.html">https://www.elabscience.com/List-detail-5594.html</a>
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? <a href="https://www.elabscience.com/List-detail-459742.html">https://www.elabscience.com/List-detail-459742.html</a>

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: <a href="www.elabscience.com">www.elabscience.com</a>
 Email: <a href="techsupport@elabscience.com">techsupport@elabscience.com</a>