Elabscience®

APC Anti-Human CD95/Fas Antibody[DX2]

Catalog No.E-AB-F1168EStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityHumanApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names	APT1,FAS1,TNFRSF6,Apo-1 antigen,FASLG receptor
Uniprot ID	P25445
Background	CD95 is a 45 kD single chain type I glycoprotein also known as Fas, APO-1, and TNFRSF6. It is
	a member of the TNF receptor superfamily. CD95 is expressed on T and B lymphocytes,
	monocytes, neutrophils, and fibroblasts. CD95 expression is upregulated by activation. The
	extracellular region of CD95 binds to CD178 (Fas ligand). CD178 binding to CD95 induces
	apoptosis and has been shown to play a role in the maintenance of peripheral tolerance.

Product Details

Form	Liquid
Size	20Tests/100Tests/100Tests×2
Clone No.	DX2
Host	Mouse
Isotype	Mouse IgG1, ĸ
Reactivity	Human
Application	FCM
Isotype Control	APC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792E]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

For Research Use Only

Elabscience®

Fluorophore


Conjugation: APC

APC is designed to be excited by the Red (627-640 nm) laser and detected using an optical filter centered near 660 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

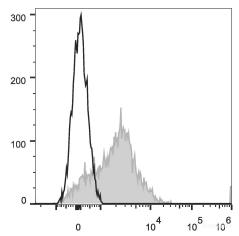

Human peripheral blood lymphocytes are stained with APC Anti-Human CD95 Antibody (filled gray histogram). Unstained lymphocytes (empty black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry <u>https://www.elabscience.com/List-detail-5594.html</u>
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>