

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

# PerCP/Cyanine5.5 Anti-Human CD11b Antibody[ICRF44]

Catalog No.E-AB-F1146JReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

### **Antigen Information**

Alternate Names Integrin alpha-M,Itgam,CD11 antigen-like family member B,CR-3 alpha chain,Leukocyte

adhesion receptor MO1,CD11b

Uniprot ID P11215

**Background** CD11b is a 165-170 kD type I transmembrane glycoprotein also known as αMintegrin, Mac-1,

CR3, and C3biR. CD11b non-covalently associates with integrin β2(CD18) and is expressed on granulocytes, monocytes/macrophages, dendritic cells, NK cells, and subsets of T and B cells. CD11b/CD18 is critical for the transendothelial migration of monocytes and neutrophils. It is also involved in granulocyte adhesion, phagocytosis, and neutrophil activation. CD11b/CD18 interacts with ICAM-1 (CD54), ICAM-2 (CD102), ICAM-4, CD14, CD23, heparin, iC3b, fibrinogen, and

factor X.

#### **Product Details**

Form Liquid

Size 20Tests/100Tests/100Tests×2

**Reactivity** Human **Application** FCM

Isotype ControlPerCP/Cyanine5.5 Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792J]Storage BufferPhosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at  $2\hbox{--}8^\circ C$  and protected from prolonged exposure to light. Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: <a href="mailto:www.elabscience.com">www.elabscience.com</a>
 Email: <a href="mailto:techsupport@elabscience.com">techsupport@elabscience.com</a>

A Reliable Research Partner in Life Science and Medicine

## **Fluorophore**


Conjugation: PerCP/Cyanine5.5

PerCP/Cyanine5.5 is designed to be excited by the blue laser (488 nm) and detected using an optical filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).

### **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

#### **Product data**

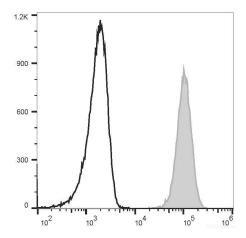

Human peripheral blood granulocytes are stained with PerCP/Cyanine5.5 Anti-Human CD11b Antibody[ICRF44] (filled gray histogram) or PerCP/Cyanine5.5 Mouse IgG1, κ Isotype Control (empty black histogram).

#### **Related Information**

- 1. Sample Preparation for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5594.html">https://www.elabscience.com/List-detail-5594.html</a>
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? <a href="https://www.elabscience.com/List-detail-459742.html">https://www.elabscience.com/List-detail-459742.html</a>

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: <a href="www.elabscience.com">www.elabscience.com</a>
 Email: <a href="techsupport@elabscience.com">techsupport@elabscience.com</a>