# **Elabscience**®

# Elab Fluor<sup>®</sup> 488 Anti-Mouse CD16/32 Antibody[2.4G2]

Catalog No.E-AB-F0997ULStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityMouseApplicationsFCM

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.

#### **Antigen Information**

Alternate Names	CD16a/b,CD32,CD32A/B,Fc fragment of IgG low affinity IIa/b receptor,Fc fragment of IgG low affinity IIIa/b receptor,Fc fragment of IgG low affinity IIIb receptor,Fc gamma receptor III A/B,FCG2A,FcGR,FCGR2A/BFCGR3,FCGR3A/B,Fc gamma RIIa/b
Uniprot ID	P08508,P08101
Background	CD16 is low affinity IgG Fc receptor III (FcR III) and CD32 is FcR II. CD16/CD32 are expressed on B cells, monocytes/macrophages, NK cells, granulocytes, mast cells, and dendritic cells. The Fc receptors bind antibody-antigen immune complexes and mediate adaptive immune responses.

#### **Product Details**

Form	Liquid
Concentration	0.5 mg/mL
Size	25µg/100µg
Clone No.	2.4G2
Host	Rat
Isotype	Rat IgG2b, κ
Reactivity	Mouse
Application	FCM
Isotype Control	Elab Fluor <sup>®</sup> 488 Rat IgG2b, κ Isotype Control[LTF-2] [Product E-AB-F09843L]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

For Research Use Only

# **Elabscience**®

## Fluorophore


#### **Conjugation:** Elab Fluor<sup>®</sup> 488

Elab Fluor<sup>®</sup> 488 is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 520 nm (e.g., a 525/40 nm bandpass filter).

#### **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the reagent to obtain optimal results [The recommended concentration is  $0.1-1 \mu g/10^6$  cells in  $100 \mu L$  volume].

### **Product data**

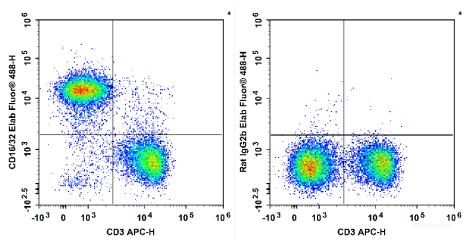

C57BL/6 murine splenocytes are stained with APC Anti-Mouse CD3 Antibody and Elab Fluor<sup>®</sup> 488 Anti-Mouse CD16/32 Antibody[2.4G2] (Left). Splenocytes are stained with APC Anti-Mouse CD3 Antibody and Elab Fluor<sup>®</sup> 488 Rat IgG2b, κ Isotype Control (Right).

### **Related Information**

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry <u>https://www.elabscience.com/List-detail-5568.html</u>
- 3. Flow Cytometry Troubleshooting Tips <u>https://www.elabscience.com/List-detail-5593.html</u>

4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>