

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 488 Anti-Mouse CD31 Antibody[390]

Catalog No.E-AB-F1180ULReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Pecam, Pecam-1,PECAM-1,CD31

Uniprot ID Q08481

Background CD31 is a 130-140 kD glycoprotein, also known as platelet endothelial cell adhesion molecule

(PECAM-1) and EndoCAM. It is a member of the Ig superfamily, expressed on endothelial cells, platelets, granulocytes, monocytes/macrophages, dendritic cells, and T and B cell subsets, and is critical for cell-cell interactions. The primary ligands for CD31 have been reported to be CD38 and the vitronectin receptor (αv $\beta 3$ integrin, CD51/CD61). Other reported functions of CD31 are

neutrophil emigration to sites of inflammation and angiogenesis.

Product Details

 Form
 Liquid

 Concentration
 0.5 mg/mL

 Size
 25μg/100μg

 Clone No.
 390

Host Rat
Isotype Rat IgG2a, κ

Reactivity Mouse
Application FCM

Isotype Control Elab Fluor[®] 488 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09833L]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® 488

Elab Fluor $^{\$}$ 488 is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 520 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use. We suggest each investigator should titrate the reagent to obtain optimal results [The recommended concentration is $0.1-1~\mu g/10^6$ cells in $100~\mu L$ volume].

Product data

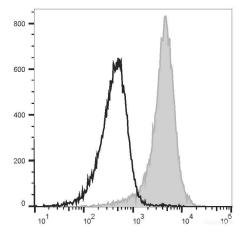

C57BL/6 murine splenocytes are stained with Elab Fluor[®] 488 Anti-Mouse CD31 Antibody[390] (filled gray histogram) or Elab Fluor[®] 488 Rat IgG2a, κ Isotype Control (empty black histogram).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086

Web: www.elabscience.com
Email: techsupport@elabscience.com