



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® 647 Anti-Human CD95/Fas Antibody[DX2]

Catalog No. E-AB-F1168M Reactivity Human Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names APT1,FAS1,TNFRSF6,Apo-1 antigen,FASLG receptor

Uniprot ID P25445

Background CD95 is a 45 kD single chain type I glycoprotein also known as Fas, APO-1, and TNFRSF6. It is

a member of the TNF receptor superfamily. CD95 is expressed on T and B lymphocytes, monocytes, neutrophils, and fibroblasts. CD95 expression is upregulated by activation. The extracellular region of CD95 binds to CD178 (Fas ligand). CD178 binding to CD95 induces apoptosis and has been shown to play a role in the maintenance of peripheral tolerance.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. DX2 Host Mouse

Mouse IgG1, κ **Isotype**

Reactivity Human **Application FCM**

Elab Fluor[®] 647 Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792M] **Isotype Control Storage Buffer** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Email: techsupport@elabscience.com Web: www.elabscience.com

Tel: 1-832-243-6086 Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® 647

Elab Fluor $^{\odot}$ 647 is designed to be excited by the Red laser (627-640 nm) and detected using an optical filter centered near 670 nm (e.g., a 660/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

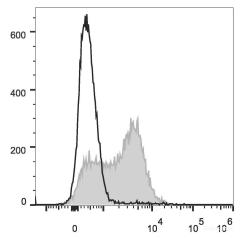

Human peripheral blood lymphocytes are stained with Elab Fluor[®] 647 Anti-Human CD95 Antibody (filled gray histogram) or Elab Fluor[®] 647 Mouse IgG1, κ Isotype Control (empty black histogram).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com