Elabscience®

Elab Fluor[®] 488 Anti-Human CD127/IL-7RA Antibody[A019D5]

Catalog No.E-AB-F1152LStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityHumanApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names	Interleukin-7 receptor subunit alpha,IL-7RA,CDw127,CD127,IL-7Ra
Uniprot ID	P16871
Background	CD127 is a 60-90 kD type I transmembrane glycoprotein also known as IL-7 receptor α chain or IL-7R α . It forms a heterodimer with the common γ chain (γ c or CD132) which is shared with the receptors for IL-2, IL-4, IL-9, IL-13, IL-15, and IL-21. CD127 is expressed on immature B cells through early pre-B stage cells, thymocytes (except CD4/CD8 double positive thymocytes), peripheral T cells, and bone marrow stromal cells. CD127 has been reported to be a useful marker for identifying memory and effector T cells. Studies have shown that CD127 expression is down-modulated on Treg cells. It can be used as a marker for differentiation of Treg and conventional T cells. The ligation of IL-7 with its receptor is important for stimulation of mature and immature T
	e e

Product Details

Form	Liquid
Size	20Tests/100Tests/100Tests×2
Clone No.	A019D5
Host	Mouse
Isotype	Mouse IgG1, ĸ
Reactivity	Human
Application	FCM
Isotype Control	Elab Fluor [®] 488 Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792L]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

For Research Use Only

Elabscience®

Fluorophore


Conjugation: Elab Fluor[®] 488

Elab Fluor[®] 488 is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 520 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

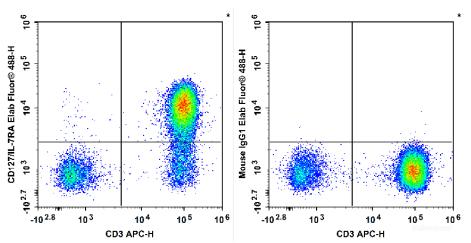

Human peripheral blood lymphocytes are stained with APC Anti-Human CD3 Antibody and Elab Fluor[®] 488 Anti-Human CD127/IL-7RA Antibody[A019D5] (Left). Lymphocytes are stained with APC Anti-Human CD3 Antibody and Elab Fluor[®] 488 Mouse IgG1, κ Isotype Control (Right).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry <u>https://www.elabscience.com/List-detail-5568.html</u>
- 3. Flow Cytometry Troubleshooting Tips <u>https://www.elabscience.com/List-detail-5593.html</u>

4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>