



A Reliable Research Partner in Life Science and Medicine

PerCP/Cyanine5.5 Anti-Mouse IgM Antibody[RMM-1]

Catalog No. E-AB-F1190J **Reactivity** Mouse **Storage** Store at 2~8°C, Avoid freeze / thaw cycles **Applications** FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Immunoglobulin heavy constant mu,IGHM,Immunoglobulin M

Uniprot ID P01872

Background IgM is the first immunoglobulin made by B cells in the immune response. Surface IgM is

expressed on the majority of mature B cells.

Product Details

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. RMM-1 Host Rat

IsotypeRat IgG2a, κReactivityMouseApplicationFCM

Isotype Control PerCP/Cyanine 5.5 Rat IgG2a, κ Isotype Control[2A3] [Product E-AB-F09832J]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com

Fax: 1-832-243-6017

Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: PerCP/Cyanine5.5

PerCP/Cyanine5.5 is designed to be excited by the blue laser (488 nm) and detected using an optical filter centered near 675 nm (e.g., a 690/50 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

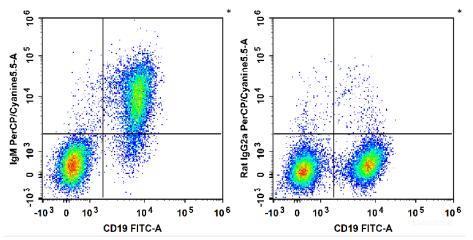

C57BL/6 murine splenocytes are stained with FITC Anti-Mouse CD19 Antibody and PerCP/Cyanine5.5 Anti-Mouse IgM Antibody (Left). Splenocytes are stained with FITC Anti-Mouse CD19 Antibody and PerCP/Cy5.5 Rat IgG2a, κ Isotype Control (Right).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com