

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

FITC Anti-Human CD122/IL-2RB Antibody[TU27]

Catalog No. E-AB-F1151C Reactivity Human Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Interleukin-2 receptor subunit beta,IL-2RB,p75,CD122,IL-2Rβ

Uniprot ID P14784

Background CD122 is a 70-75 kD type I transmembrane glycoprotein and member of the Ig superfamily. It is

IL-2 receptor β chain also known as IL-2Rβ, which is also shared by the IL-15 receptor. CD122 is constitutively expressed by NK cells and at lower levels by a subset of T cells. Its expression is upregulated upon activation. The IL-2R β chain can combine with either the common γ subunit (yc, CD132) alone or with the yc subunit and the IL- $2R\alpha$ subunit (CD25) to generate intermediate

or high affinity IL-2 receptor complexes, respectively. CD122 expression levels can be

upregulated by activation.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. **TU27** Host Mouse

Mouse IgG1, κ **Isotype**

Reactivity Human **Application FCM**

Isotype Control FITC Mouse IgG1, K Isotype Control[MOPC-21] [Product E-AB-F09792C]

Storage Buffer Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C Stability & Storage Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Email: techsupport@elabscience.com Web: www.elabscience.com

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

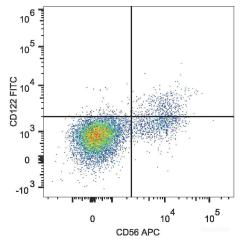

Human peripheral blood lymphocytes are stained with FITC Anti-Human CD122/IL-2RB Antibody and APC Anti-Human CD56 Antibody.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com