

Fax: 1-832-243-6017



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® Violet 450 Anti-Mouse CD11c Antibody[N418]

E-AB-F0991Q Catalog No. Reactivity Mouse Storage Store at 2~8°C, Avoid freeze / thaw cycles **Applications FCM**

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names Integrin alpha-X,Itgax,CD11 antigen-like family member C,Leukocyte adhesion receptor

p150+95,CD11c

Uniprot ID Q9QXH4

Background CD11c is a 150 kD glycoprotein also known as αX integrin, CR4, and p150. CD11c forms a

> αΧβ2 heterodimer with β2 integrin (CD18). It is primarily expressed on dendritic cells, NK cells, a subset of intestinal intraepithelial lymphocytes (IEL), and some activated T cells. The $\alpha X\beta 2$ integrin plays an important role in cell-cell contact by binding its ligands: iC3b, fibrinogen and

CD54.

Product Details

Form Liquid

50Tests/100Tests/100Tests×2 Size

Clone No. N418

Armenian Hamster Host **Isotype** Armenian Hamster IgG

Mouse Reactivity **FCM Application**

Elab Fluor[®] Violet 450 Armenian Hamster IgG Isotype Control[PIP] [Product E-AB-F09852Q] **Isotype Control** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. **Storage Buffer**

Biological ice pack at 4 °C Shipping Keep as concentrated solution. Stability & Storage

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Email: techsupport@elabscience.com Web: www.elabscience.com

Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® Violet 450

Elab Fluor[®] Violet 450 is designed to be excited by the violet laser (405 nm) and detected using an optical filter centered near 450 nm (e.g., a 450/45 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

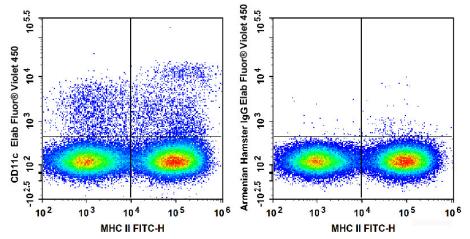

C57BL/6 murine splenocytes are stained with FITC Anti-Mouse MHC II Antibody and Elab Fluor[®] Violet 450 Anti-Mouse CD11c Antibody (Left). Splenocytes are stained with FITC Anti-Mouse MHC II Antibody and Elab Fluor[®] Violet 450 Armenian Hamster IgG Isotype Control (Right).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com