



A Reliable Research Partner in Life Science and Medicine

Elab Fluor® Violet 450 Anti-Mouse CD4 Antibody[GK1.5]

Catalog No.E-AB-F1097QReactivityMouseStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names T-cell surface glycoprotein CD4,CD4,T-cell surface antigen T4/Leu-3,CD4

Uniprot ID P06332 **Gene ID** 12504

Background CD4 is a 55 kD protein also known as L3T4 or T4. It is a member of the Ig superfamily,

primarily expressed on most thymocytes, a subset of T cells, and weakly on macrophages and dendritic cells. It acts as a coreceptor with the TCR during T cell activation and thymic differentiation by binding MHC class II and associating with the protein tyrosin kinase, lck.

Product Details

Form Liquid

Size 50Tests/100Tests/100Tests×2

Clone No. GK1.5 Host Rat

IsotypeRat IgG2b, κReactivityMouseApplicationFCM

Isotype ControlElab Fluor® Violet 450 Rat IgG2b, κ Isotype Control[LTF-2] [Product E-AB-F09842Q]Storage BufferPhosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086
Web: www.elabscience.com Email: techsupport@elabscience.com

Tel: 1-832-243-6086 Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: Elab Fluor® Violet 450

Elab Fluor[®] Violet 450 is designed to be excited by the violet laser (405 nm) and detected using an optical filter centered near 450 nm (e.g., a 450/45 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

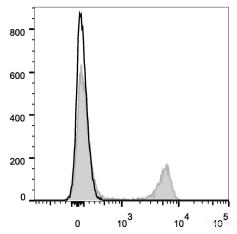

Mouse splenocytes are stained with Elab Fluor[®] Violet 450 Anti-Mouse CD4 Antibody (filled gray histogram). Unstained splenocytes (blank black histogram) are used as control.

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com