

#### A Reliable Research Partner in Life Science and Medicine

# Recombinant Mouse FGCP/GCPII/mGCP/NAALADase I/PSMA(N-6His)

Catalog No. PKSM041421

Note: Centrifuge before opening to ensure complete recovery of vial contents.

## **Description**

Synonyms Glutamate carboxypeptidase 2;FGCP;GCPII;mGCP;NAALADase I;PSMA;Cell

growth-inhibiting gene 27 protein; Folate hydrolase 1

**Species** Mouse

Expression Host HEK293 Cells
Sequence Ile44-Ala752
Accession Q0VDM5
Calculated Molecular Weight 80.7 kDa
Observed molecular weight 100-120 kDa
Tag N-His

**Bioactivity** Not validated for activity

## **Properties**

**Purity** > 95 % as determined by reducing SDS-PAGE.

**Endotoxin** < 1.0 EU per µg of the protein as determined by the LAL method.

**Storage** Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to

-80°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots

of reconstituted samples are stable at < -20°C for 3 months.

**Shipping** This product is provided as lyophilized powder which is shipped with ice packs.

Formulation Lyophilized from a 0.2 µm filtered solution of 20mM PB, 8% Trehalose, 4%

Mannitol, 50mM NaCl, 0.05% Tween 80, pH7.5.


Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as

protectants before lyophilization.

Please refer to the specific buffer information in the printed manual.

**Reconstitution** Please refer to the printed manual for detailed information.

## Data

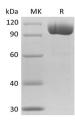

> 95 % as determined by reducing SDS-PAGE.

# Background

### For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017

Web: <u>www.elabscience.com</u> Email: <u>techsupport@elabscience.com</u>

## **Elabscience Bionovation Inc.**



A Reliable Research Partner in Life Science and Medicine

Glutamate carboxypeptidase 2, also known as FOLH1, PSMA, belongs to the M28B subfamily and the peptidase M28 family. It is highly expressed in prostate epithelium and can be detected in urinary bladder, kidney, testis, ovary, fallopian tube, breast, adrenal gland, liver, esophagus, stomach, small intestine, colon and brain (at protein level). PSMA is used as a diagnostic and prognostic indicator of prostate cancer, and as a possible marker for various neurological disorders such as schizophrenia, Alzheimer disease and Huntington disease. It has both folate hydrolase and N-acetylated-alpha-linkedacidic dipeptidase (NAALADase) activity and has a preference for tri-alpha-glutamate peptides. PSMA involves in prostate tumor progression and also exhibits a dipeptidyl-peptidase IV type activity. In vitro, PSMA cleaves Gly-Pro-AMC. PSMA is stable at pH greater than 6.5.

For Research Use Only

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017 Email: techsupport@elabscience.com

Web: www.elabscience.com