



A Reliable Research Partner in Life Science and Medicine

PE/Elab Fluor® 594 Anti-Human CD44 Antibody[P2A1]

Catalog No.E-AB-F1038PReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names CD44 antigen, CD44, CDw44, Epican, Phagocytic glycoprotein 1, PGP-1, Phagocytic glycoprotein

I,PGP-I,CD44,LHR, MDU2, MDU3, MIC4

Uniprot ID P16070 Gene ID 960

Background CD44 is a 80-95 kD glycoprotein also known as Hermes, Pgp1, H-CAM, or HUTCH. It is

expressed on all leukocytes, endothelial cells, hepatocytes, and mesenchymal cells. As B and T cells become activated or progress to the memory stage, CD44 expression increases from a low or mid level of intensity to high expression levels. Thus, CD44 has been reported to be a valuable marker for memory cell subsets. CD44 is an adhesion molecule involved in leukocyte attachment to and rolling on endothelial cells, homing to peripheral lymphoid organs and to the sites of

inflammation, and leukocyte aggregation.

Product Details

Form Liquid

Size 20Tests/100Tests/100Tests×2

Clone No. P2A1 Host Mouse

Isotype Mouse IgG2a, κ

Reactivity Human **Application** FCM

Isotype Control Storage BufferPE/Elab Fluor[®] 594 Mouse IgG2a, κ Isotype Control[C1.18.4] [Product E-AB-F09802P]

Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

Shipping Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

Web: www.elabscience.com Email: techsupport@elabscience.com

Toll-free: 1-888-852-8623 Tel: 1-832-243-6086 Fax: 1-832-243-6017

A Reliable Research Partner in Life Science and Medicine

Fluorophore


Conjugation: PE/Elab Fluor® 594

PE/Elab Fluor[®] 594 is designed to be excited by the blue (488 nm), Green (532 nm) and yellow-green (561 nm) lasers and detected using an optical filter centered near 620 nm (e.g., a 610/20 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

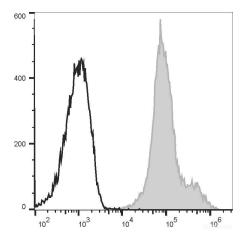

Human peripheral blood lymphocytes are stained with PE/Elab Fluor[®] 594 Anti-Human CD44 Antibody[P2A1] (filled gray histogram) or PE/Elab Fluor[®] 594 Mouse IgG2a, κ Isotype Control (empty black histogram).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? https://www.elabscience.com/List-detail-459742.html

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: www.elabscience.com
 Email: techsupport@elabscience.com