Elabscience®

FITC Anti-Human CD185/CXCR5 Antibody[J252D4]

Catalog No.E-AB-F1287CStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityHumanApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

| Alternate Names | CXCR5,C-X-C chemokine receptor type 5,Burkitt lymphoma receptor 1,BLR1,Monocyte- |
|-----------------|---|
| | derived receptor 15,MDR15 |
| Uniprot ID | P32302 |
| Gene ID | 643 |
| Background | CD185, also known as CXCR5, is a 42 kD G-protein coupled receptor with seven transmembrane regions. CXCR5 is expressed by mature B cells, follicular helper T cells, Burkitt's lymphoma cells and a subset of neurons, and mediates cell migration to the B cell follicles in the secondary lymphoid organs. The ligand of CXCR5 is CXCL13 (BLC). |

Product Details

| Form | Liquid |
|---------------------|---|
| Size | 20Tests/100Tests/100Tests×2 |
| Clone No. | J252D4 |
| Host | Mouse |
| Isotype | Mouse IgG1, ĸ |
| Reactivity | Human |
| Application | FCM |
| Isotype Control | FITC Mouse IgG1, κ Isotype Control[MOPC-21] [Product E-AB-F09792C] |
| Storage Buffer | Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant. |
| Shipping | Biological ice pack at 4 °C |
| Stability & Storage | Keep as concentrated solution. |
| | Store at 2~8°C and protected from prolonged exposure to light.Do not freeze. |
| | This product is guaranteed up to one year from purchase. |

For Research Use Only

Elabscience®

Fluorophore


Conjugation: FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

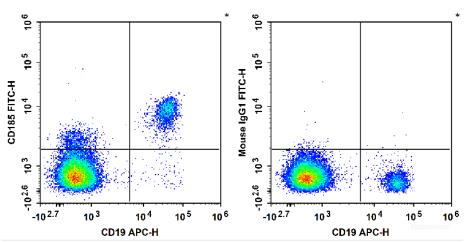

Human peripheral blood lymphocytes are stained with APC Anti-Human CD19 Antibody and FITC Anti-Human CD185/CXCR5 Antibody (Left). Lymphocytes are stained with APC Anti-Human CD19 Antibody and FITC Mouse IgG1, κ Isotype Control (Right).

Related Information

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry <u>https://www.elabscience.com/List-detail-5568.html</u>
- 3. Flow Cytometry Troubleshooting Tips <u>https://www.elabscience.com/List-detail-5593.html</u>

4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>