Elabscience®

PE Anti-Human CD49d Antibody[HP1/2]

Catalog No.E-AB-F1359DStorageStore at 2~8°C, Avoid freeze / thaw cycles

Reactivity Applications

Human FCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

Antigen Information

Alternate Names	VLA-4 α chain, α 4 integrin, Integrin α 4 chain, ITGA4
Uniprot ID	P13612
Gene ID	3676
Background	CD49d is a 150 kD α integrin chain known as α 4 integrin or VLA-4 α chain. It forms a
	heterodimer with either integrin β 1 (α 4 β 1, VLA-4) or β 7 (α 4 β 7). CD49d is expressed broadly on T lymphocytes, B lymphocytes, monocytes, thymocytes, eosinophils, basophils, mast cells, NK cells, dendritic cells, and some non-hematopoietic cells, but not on normal red blood cells, platelets or neutrophils. VLA-4 binds to VCAM-1 (CD106) and fibronectin. α 4 β 7 is the receptor for VCAM-1 and MAdCAM-1. CD49d participates in mononuclear cell trafficking to endothelial sites of inflammation and has roles in cell-cell interactions and cell adhesion to extracellular matrices. CD49d is involved in lymphocyte migration, T cell activation, and hematopoietic stem cell differentiation. CD49d is a marker to isolate pure populations of Treg cells due to its absence on Foxp3+ cells.

Product Details

Form	Liquid
Size	20Tests/100Tests/100Tests×2
Clone No.	HP1/2
Host	Mouse
Isotype	Mouse IgG2a, ĸ
Reactivity	Human
Application	FCM
Isotype Control	PE Mouse IgG2a, κ Isotype Control[C1.18.4] [Product E-AB-F09802D]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

For Research Use Only

Elabscience®

Fluorophore


Conjugation: PE

PE is designed to be excited by the Blue (488 nm), Green (532 nm) and Yellow-Green (561 nm) lasers and detected using an optical filter centered near 575 nm (e.g., a 585/42 nm bandpass filter).

Recommended usage

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5 μ L of antibody per test (million cells in 100 μ L staining volume or per 100 μ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

Product data

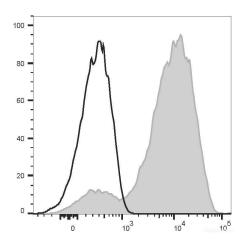

Human peripheral blood lymphocytes are stained with PE Anti-Human CD49d Recombinant Antibody (filled gray histogram) or PE Mouse IgG2a, κ Isotype Control (empty black histogram).

Related Information

- 1. Sample Preparation for Flow Cytometry <u>https://www.elabscience.com/List-detail-5594.html</u>
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips https://www.elabscience.com/List-detail-5593.html
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>