



A Reliable Research Partner in Life Science and Medicine

## FITC Anti-Human CD123 Antibody[HI12H7]

Catalog No.E-AB-F1340CReactivityHumanStorageStore at 2~8°C, Avoid freeze / thaw cyclesApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

### **Antigen Information**

**Alternate Names** IL-3Rα, IL-3 Receptor alpha

**Uniprot ID** P26951 **Gene ID** 3563

**Background** CD123 is the 70 kD transmembrane α chain of the IL-3 receptor. Alone, CD123 binds IL-3 with

low affinity; when CD123 associates with CDw131 (common  $\beta$  chain), it binds IL-3 with high affinity. CD123 does not transduce intracellular signals upon binding IL-3 and requires the  $\beta$  chain for this function. CD123 is expressed by myeloid precursors, macrophages, dendritic cells,

mast cells, basophils, megakaryocytes, and some B cells.

#### **Product Details**

Form Liquid

Size 20Tests/100Tests×2

Clone No. HI12H7 Host Mouse

**Isotype** Mouse IgG1,  $\kappa$ 

**Reactivity** Human **Application** FCM

**Isotype Control** FITC Mouse IgG1,  $\kappa$  Isotype Control[MOPC-21] [Product E-AB-F09792C]

**Storage Buffer** Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.

**Shipping** Biological ice pack at 4 °C **Stability & Storage** Keep as concentrated solution.

Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.

This product is guaranteed up to one year from purchase.

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086
 Fax: 1-832-243-6017

 Web: <a href="https://www.elabscience.com">www.elabscience.com</a>
 Email: <a href="mailto:techsupport@elabscience.com">techsupport@elabscience.com</a>

A Reliable Research Partner in Life Science and Medicine

# **Fluorophore**


**Conjugation: FITC** 

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

# **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

### **Product data**

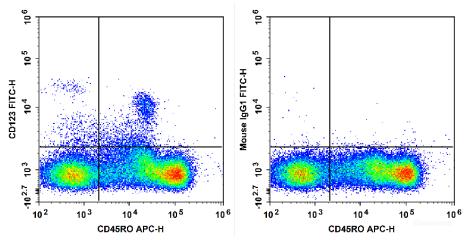

Human peripheral blood mononuclear cells are stained with APC Anti-Human CD45RO Antibody and FITC Anti-Human CD123 Antibody[HI12H7] (Left). Mononuclear cells are stained with APC Anti-Human CD45RO Antibody and FITC Mouse IgG1,  $\kappa$  Isotype Control (Right).

### **Related Information**

- 1. Sample Preparation for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5594.html">https://www.elabscience.com/List-detail-5594.html</a>
- 2. Staining Cell Surface Targets for Flow Cytometry <a href="https://www.elabscience.com/List-detail-5568.html">https://www.elabscience.com/List-detail-5568.html</a>
- 3. Flow Cytometry Troubleshooting Tips <a href="https://www.elabscience.com/List-detail-5593.html">https://www.elabscience.com/List-detail-5593.html</a>
- 4. How to select the appropriate detection channel through the spectrogram? <a href="https://www.elabscience.com/List-detail-459742.html">https://www.elabscience.com/List-detail-459742.html</a>

For Research Use Only

 Toll-free: 1-888-852-8623
 Tel: 1-832-243-6086

 Web: <a href="www.elabscience.com">www.elabscience.com</a>
 Email: <a href="techsupport@elabscience.com">techsupport@elabscience.com</a>